visomat

double comfort

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Upper Arm Blood Pressure Monitor

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Table of contents

A	Intended use	38
B	Safety instructions	40
C	Operating the device	
1.	Device description	43
2.	Display	44
3.	Important usage information	45
4.	Initial operation of the device	45
5.	Inserting/replacing batteries, battery icon	46
6.	Activating and setting the date/time	46
7.	Attaching the cuff	48
8.	Body position during the measurement	49
9.	Selecting the user	50
10.	Measuring blood pressure	50
11.	Classification of the measured values (WHO)	53
12.	Pulse pressure	54
13.	Cardiac arrhythmia/atrial fibrillation/pulse disturbances	54
14.	Motion control	55
15.	Using the memory	55
16.	Transferring the data to a PC	57
D	What you should know about blood pressure	
1.	The systolic and diastolic blood pressure value	58
2.	Why you measure different values	58
3.	Why should you measure blood pressure regularly?	58
E	Technical information	
1.	Error messages	59
2.	Customer service	61
3.	Technical data	61
4.	Original spare parts and accessories	62
5.	Metrological inspection	63
6.	Explanation of symbols	63
7.	Disposal	65

Table of contents

F	Maintaining the device	66
G	Warranty	67

Thank you for choosing the visomat double comfort upper arm blood pressure monitor.

These instructions for use are intended to help the user operate the device safely and efficiently. The unit must be used in accordance with the procedures contained in these instructions for use and must not be used for any other purposes.

It is important that you read all the instructions carefully before using the device. Please pay particular attention to the chapter "Important usage information" on page 45.

This device is intended for laypersons aged 12 or over for performing non-invasive self-measurement of systolic and diastolic blood pressure, determining pulse rate and calculating pulse pressure on the upper arm.

How it works

The device determines the blood-pressure values simultaneously using two different measurement procedures, the Korotkoff procedure and the oscillometric measurement procedure. As a result, the measurement values determined are very reliable and precise.

In the Korotkoff procedure, a highly sensitive microphone in the cuff is used to record the sounds that arise during the blood pressure measurement due to the turbulence of the blood in the artery. The stethoscope blood pressure measurement performed by doctors also functions according to this principle.

The oscillometric measurement procedure that is customary in fully automatic blood pressure measurement devices uses a sensor to record the small fluctuations in pressure (oscillations) inside the cuff that are caused by the expansion and contraction of the arteries in the arm (heartbeat).

visomat double comfort uses both procedures to determine the measurement values and indicates the result on the display.

Indications

Hypertension, hypotension, tachycardia, bradycardia, irregular pulses.

Contraindications

- During pregnancy and in cases of pre-eclampsia, the device may only be used after consulting your doctor.
- Do not under any circumstances place the cuff over a critical area, e.g. a wound, aneurysm, etc. or on an arm with an arteriovenous shunt. Risk of injury! Any supply via an intravascular access point (infusion) or other medical monitoring devices could possibly be interrupted.

- Do not use the device without first consulting your doctor if you are having dialysis treatment or are taking anticoagulants, platelet aggregation inhibitors or steroids. These circumstances may cause internal bleeding.
- Do not take a blood pressure measurement on the same side as a mastectomy with simultaneous removal of the axilla lymph nodes.
- This device must not be used by persons who are unable to follow the instructions for use, operate the device properly or read the measurement results accurately.
- Do not use the device on children.

Safety instructions


- The cuff must only be used on the upper arm and must never be used on other limbs.
- The device is not a diagnostic instrument. It only indicates the potential presence of blood pressure-related conditions. Discuss the measurements with your doctor to obtain an accurate diagnosis.
- The display of the pulse frequency is not suitable for checking the frequency of cardiac pacemakers. Cardiac pacemakers and blood pressure monitors do not influence each other in terms of their mode of operation.
- The device is not suitable for use in close proximity to electrosurgical instruments or magnetic resonance imaging equipment.
- Do not use the device in explosive environments such as in the vicinity of flammable anaesthetics or in an oxygen chamber.
- Ensure that the air hose is not kinked. A kinked air hose may prevent the cuff from deflating, thus interrupting the blood flow in the arm for too long.
- Please wait a few minutes between successive measurements, otherwise the blood flow in the arm is interrupted for too long and injuries could occur. For this purpose, use the triple measurement function, which automatically controls and alternates the measurements and the intervals for recovering blood circulation (see page 52).
- The device contains small parts that could be swallowed by children. The air hose poses a strangulation hazard. Do not leave the device unattended with children.
- Discuss the measured values with your doctor. Do not assess the measurement results yourself. Never change the dose of medicines prescribed by your doctor.
- Before carrying out your self-measurement, please pay attention to the chapter "Important usage information" on page 45.

Safety instructions

- The unit contains sensitive parts and must be protected from extreme temperature fluctuations, humidity, shocks, dust and direct sunlight.
- If the device has been stored at temperatures below 10 °C or above 40 °C, allow it to stand at room temperature for at least 2 hours before using it.
- Do not knock or drop the device.
- The device is not waterproof. Never immerse it in water or any other liquids. If liquid enters the device, this may result in serious damage and may cause malfunctions.
- Do not apply pressure to the display. Do not place the device with the display facing downwards.
- The unit must only be operated with the cuff and accessories provided (see "Original spare parts and accessories" on page 62). Using different cuffs or accessories can lead to incorrect measurement results. The warranty will be invalidated if the unit is damaged by third-party accessories!
- A consistently good power supply to the device is necessary for fault-free blood pressure measurements.
 - Only use long-life alkaline batteries (LR6).
 - Always replace all batteries simultaneously during battery replacement. Do not mix new and old batteries or batteries of different types.
 - You will require 4 x 1.5 V batteries. Rechargeable batteries only have a voltage of 1.2 V and are therefore not suitable.
 - When using the device with a mains adapter, please only use the visomat mains adapter, which has been specially tested for medical devices.
- Never simultaneously touch the device's live electrical parts and the patient.


- If the device is not to be operated for an extended period of time, please remove the batteries. In principle, all batteries can leak.
- Do not remove batteries from the device or separate it from the mains adapter when the device is switched on.
- Using this device in the vicinity of mobile phones, microwaves or other devices with strong electromagnetic fields can cause malfunctions and inaccurate measurements. During use, maintain a distance of at least 30 cm from such devices.
- Never open or modify the device or the cuff (except for replacing the batteries). If the unit has been opened, it must be subjected to a metrological inspection by an authorised institution.
- Ensure that the cuff connector's microphone plug is plugged firmly into the device before you begin measuring.
- To avoid inaccurate readings, please adhere to the intended operating and storage conditions. See "Technical data" on page 61.
- The inflating and measuring procedure can be interrupted by pressing the Start/Stop button or by removing the cuff. The device then stops the inflation procedure and deflates the cuff.
- If the device malfunctions or is faulty, please refer to the troubleshooting guide starting on page 59, or contact customer service (see page 61).
- Any serious adverse effects linked to the product must be reported to the manufacturer and the relevant authorities in the member state in which the user/patient resides.

1. Device description

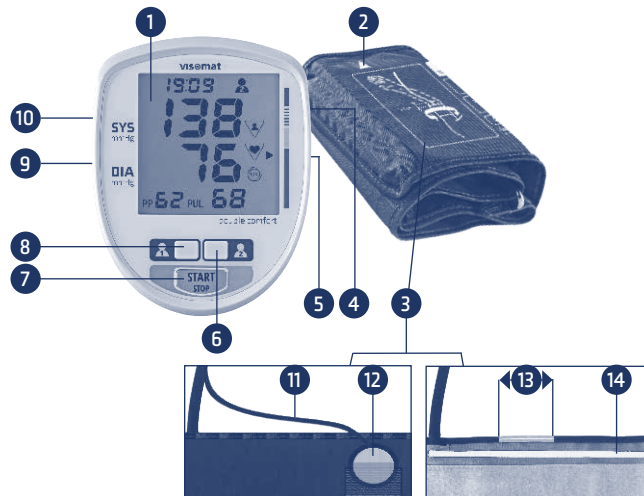- | | |
|---|---|
| 1 Display | 8 Memory button user 1 (memory recall button) |
| 2 Measurement arrow for arm circumference | 9 Battery compartment |
| 3 Cuff | 10 Cuff connection |
| 4 USB connecting socket | 11 Microphone cable |
| 5 Connecting socket for mains adapter | 12 Microphone |
| 6 Memory button user 2 (memory recall button) | 13 Artery marking |
| 7 Start/Stop button | 14 Marking for arm circumference |

2. Display

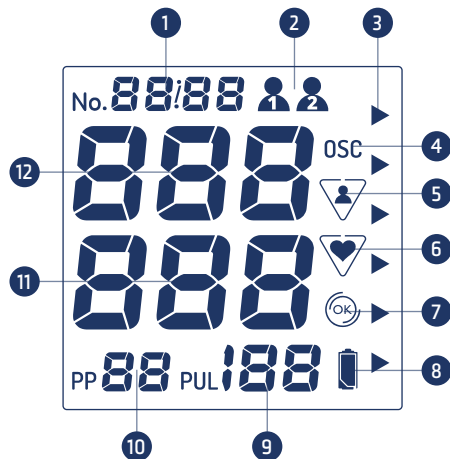- | | | | |
|---|---|----|---|
| 1 | Memory capacity or date/
time or countdown for triple
measurement | 7 | Cuff position check |
| 2 | User | 8 | Battery indicator |
| 3 | WHO classification (page 53) | 9 | PUL 1/min = pulse
Calculated pulse frequency
per minute |
| 4 | Oscillometric measurement
(page 52) | 10 | Pulse pressure |
| 5 | Motion control | 11 | DIA = diastole |
| 6 | Pulse signal display or irregu-
lar pulse rates | 12 | SYS = systole |

3. Important usage information

The measurement results of automated blood pressure monitors can be influenced by the measuring location, posture, previous exertion and general physical condition. Please pay attention to the usage information in order to obtain accurate readings.

- Do not consume alcohol, caffeine or smoke for at least one hour before measuring.
- Rest for at least 5 minutes before taking the measurement. Depending on the previous degree of exertion, you may even need to rest for up to an hour.
- Expose the upper arm. Under no circumstances may clothes hinder the blood flow to or from the arm, as this affects the blood pressure at the measuring point and can lead to incorrect readings.
- Do not move or talk during measuring.
- Breath slowly and deeply. Do not hold your breath.
- In order to exclude side differences and to obtain comparable measurement results, it is important to always take the measurement on the same arm. Ask your doctor which side it is advisable for you to measure.
- Always take your daily blood pressure measurements at around the same time.
- Blood pressure is not a fixed value. It can go up or down by more than 20 mmHg in patients within a few minutes.

4. Initial operation of the device

Insert the batteries provided into the device.

If the device is to be operated using mains power, the cable connector of the mains adapter (not supplied) must be plugged into the socket on the right-hand side of the device. The batteries are automatically switched off.

Please only use the visomat mains adapter. See “Original spare parts and accessories” on page 62.

5. Inserting/replacing batteries, battery icon

- Remove the battery compartment cover on the underside of the device.
- Remove the old batteries from the device and insert the new ones. Please ensure correct polarity (marking in battery compartment).
- Close the battery compartment by clipping the battery cover back into the device.
- When the “empty battery ” symbol first appears in the display, you are still able to take around 30 measurements. Please replace the batteries within this time frame.

6. Activating and setting the date/time

The device has a date/time function that can be activated as required. This function is deactivated on delivery.

When the date/time function is activated, the measurement results are saved along with their date and time. This can make it easier to view the readings at a later point.

If you are using the visomat double comfort software to analyse measured values on your computer, the date/time function must be activated so that the date and time can be stored with the measurements. Measured values that are stored without the date and time cannot be transferred to the computer.

Activating and deactivating the date/time function

Press the Start/Stop button to activate the date/time function. Then immediately press one of the two Memory buttons while the full display is still visible. The function can be deactivated at any time using the same procedure.

When the date/time function is activated, the device continues to display the time even when it is switched off. The power consumption for this is extremely low.

Setting the date/time

Initially, the display will show the four-digit year number (Figure 1). Press Memory button 1 to reduce the year number or Memory button 2 to increase the year number. Confirm the entry by pressing the Start/Stop button.

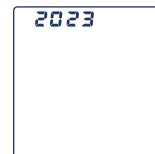
Figure 1

You are now able to set the month. (Figure 2). Press Memory button 1 to reduce the value or Memory button 2 to increase the value. Set the correct month and confirm the entry by pressing the Start/Stop button.

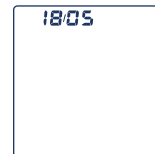
Figure 2

Repeat these steps to enter the day (Figure 2) and the hours and the minutes (Figure 3).

After replacing the batteries, you will have to reset the date and time.

Changing the date/time

With the device switched off, hold down both memory buttons simultaneously for three seconds to enable the date and time to be set. Perform the setting as described on page 46.

If you take the opportunity to evaluate your blood pressure values on a PC, the date and time are automatically compared with your PC's system time

Figure 3

during the data transfer. In case of deviations, a query will appear on the PC screen asking whether the device's date and time should be adjusted.

7. Attaching the cuff

Before applying the cuff, please check that the circumference of your upper arm is within the range specified on the cuff. An incorrect cuff size could result in inaccurate readings.

- Expose your upper arm.
- Push the cuff onto the arm until the lower edge of the cuff is 2-3 cm above the crook of the arm (Figure 1, Figure 2).
- When placed on the left arm, the air hose runs down the middle of the crook of the arm to the device. This automatically positions the 4-cm long artery marking and the microphone above the pulse position (Figure 3).
- When placed on the right arm, the cuff must be rotated to the left until the artery marking and the microphone are above the pulse position. The air hose then runs along the inner side of the upper arm (Figure 4).
- Now pull the free end of the cuff tight and close the Velcro fastener (Figure 5).
- Check that the measurement arrow on the cuff is within the "Marking for arm circumference" on the cuff edge (Figure 6).

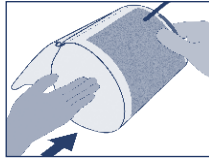
Figure 1

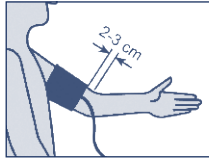
Figure 2

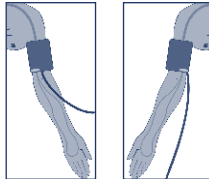
Figure 3

Figure 4

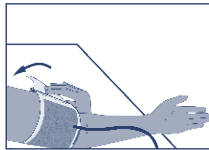
Figure 5

Figure 6

- Connect the plug of the cuff to the cuff connection on the left-hand side of the device (Figure 7). Ensure that the cuff connector's microphone plug is plugged firmly into the device.
- Place the arm with the cuff on a table so that it is loosely stretched out and hold it very still during the measurement, do not talk. The palm of your hand should be facing upwards.


Figure 7

8. Body position during the measurement

The measurement should preferably be performed in a sitting position. Measurements in a lying position are also permitted in exceptional cases. The body must be in a relaxed position.

Measurement in a sitting position

- Sit at a table (preferably at the height of a dining table, not a coffee table).
- Sit with your back against the backrest of the chair.
- Rest your entire lower arm on the table with the palm facing upwards.
- Put your feet on the floor and do not cross your legs.

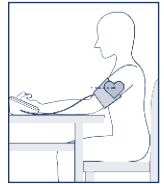Measurement in a lying position

- Lie on your back. Look up at the ceiling during the measurement. Keep still and do not move.
- Ensure that the measuring point on the upper arm is approximately level with your heart. If the measuring point is lower than the heart, the measured values may be higher. If the measuring point is higher than the heart, the measured values may be lower.

The results for measurements taken in a sitting and lying position may differ slightly.

9. Selecting the user

The device is designed for two users. The values for the two users are stored separately and can also be evaluated separately.

You may also use the two user memories to evaluate measurements taken in morning and those taken in the evening separately. To do this, carry out the measurement in the morning using user 1 and in the evening using user 2.

- Before taking the measurement, check that the correct user is displayed.
- To change the user, press the user button for whose memory bank you want to switch over to.
- Now carry out the measurement.
- Also, you will be given another opportunity to assign the measurement to the correct user after the measurement has been taken (see page 51).

10. Measuring blood pressure

You can choose to perform either one single measurement or one triple measurement, in which three automatically controlled single measurements are performed successively with brief intervals between them.

To perform a single measurement, press the start/stop button briefly. To start a triple measurement, press and hold the start/stop button for at least two seconds.

Performing a single measurement

- Press the start/stop button once briefly to switch on the device. The full display is shown (Figure 1):
- After the device has completed its calibration against the ambient air pressure (Figure 2), the automatically controlled inflation process starts. The cuff is inflated to the pressure required for measurement.
- Shortly after the inflation procedure begins, a check is carried out to ensure that the cuff has not been applied too loosely. If the cuff is applied correctly, the ☺ symbol is displayed (Figure 3). If the symbol is not displayed by the time the pulse signal starts, the cuff is too loose. In this case, you should repeat the measurement after you have applied the cuff correctly (see page 48).
- The actual measurement process then begins as the cuff deflates. The ♥ symbol flashes in the display (Figure 4). In addition, beeps denote the pulse.
- A long beep denotes the end of measurement. The cuff is automatically deflated.
- The determined systole, diastole and pulse values, as well as the pulse pressure (PP) are shown in the display (Figure 5).

When the measurement is finished, you will be given another opportunity to change the memory to which the measured value is to be saved. To do this, press the Memory button of the user for whom the measured value is to be saved.

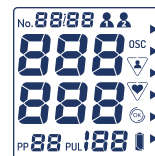
Figure 1

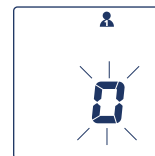
Figure 2

Figure 3

Figure 4

Figure 5

The device switches itself off after approximately 3 minutes. You can also switch the device off before this by pressing the Start/Stop button.

If the acoustic signal is too weak, the microphone integrated in the cuff will not be able to record the Korotkoff sounds of the blood flow without errors.

In this case, the measurement values are determined according to the oscillometric procedure (see page 38). The measurement values are indicated together with the symbol **OSC** (Figure 6). In this case, check whether the cuff and the microphone are correctly seated and whether the microphone jack has been fully inserted into the device. Repeat the measurement. During the measurement, the palm of the hand must be turned upwards.

Figure 6

If the symbol **OSC** appears constantly during each measurement, please contact visomat customer service (see page 61).

Performing a triple measurement

- Press and hold the start/stop button for at least two seconds to switch on the device and to start the triple measurement. The full display is shown (Figure 7). If you press and hold the start/stop button for more than 10 seconds, the measurement is stopped and the device switches off.
- The first single measurement will then begin. The continued sequence of measurement is described under "Performing a single measurement" on page 51. In addition, No. 1/3 (measurement 1 of 3) is indicated on the top left of the display (Figure 8).
- At the end of the measurement, the result is shown in the display. There is an interval of two minutes to prevent the blood flow in the arm from being interrupted for too long. A clock counts down from 120 seconds to 0 at the top of the display (Figure 9).

Figure 7

Figure 8

Figure 9

- After the 120 seconds elapses, the second measurement starts automatically. The display shows No. 2/3 (measurement 2 of 3). After the second measurement is complete, there is another interval of 120 seconds (Figure 9). Then the third measurement begins.
- Once the third measurement is complete, the average values are calculated and shown on the display (Figure 10). Only measurements 2 and 3 are included in the calculation of the average values. The first measurement is discarded and not considered in the calculation.

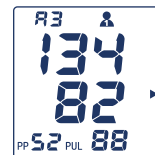
Figure 10

The device switches itself off after approximately 3 minutes. You can also switch the device off before this by pressing the Start/Stop button.

11. Classification of the measured values (WHO)

The device categorises the blood pressure readings according to the 1999 guidelines of the World Health Organisation (WHO). You can determine the classification after each measurement using the triangle in the display and the colour fields next to the display.

Classification	Systolic pressure = top value mmHg	Diastolic pressure = bottom value mmHg
Stage 3 hypertension	≥ 180	≥ 110
Stage 2 hypertension	160-179	100-109
Stage 1 hypertension	140-159	90-99
High-normal	130-139	85-89
Normal	120-129	80-84
Optimal	< 120	< 80

The WHO classification is stored along with the measured values and both can be recalled from the memory.

Depending on age, weight and general health, the blood pressure values can be different. Only a doctor can determine the right blood pressure range for you and assess whether your blood pressure has reached a dangerous level for you. Discuss your blood pressure values with your doctor. **Never change the dose of medicines prescribed by your doctor!**

12. Pulse pressure

The pulse pressure – not to be confused with the pulse beat – provides an indication of the elasticity of the blood vessels. A stiff system of vessels can have a negative effect on the cardiovascular system. Studies show that the cardiovascular risk increases if the pulse pressure value is permanently higher than 65 mmHg.

The heart works in two phases, the contraction phase (expulsion phase, systole) and the relaxation phase (filling or resting phase, diastole). The pressure difference between systole and diastole is called the pulse pressure or pulse amplitude. The higher the pulse pressure, the stiffer your vessels are thought to be.

High pulse pressure	Above 65 mmHg
Increased pulse pressure	55 to 65 mmHg
Normal pulse pressure	Below 55 mmHg

If your pulse pressure is permanently above 55 mmHg, consult your doctor.

13. Cardiac arrhythmia/atrial fibrillation/pulse disturbances

If the device detects irregular pulse rates during the measurement, the symbol is displayed once the measurement is complete. This can be triggered by cardiac arrhythmias, movement, talking or even by deep breathing.

The symbol is stored with the respective measurement.

Symbol is flashing: The device has detected strong pulse disturbances with a very fast pulse. This may be atrial fibrillation. Repeat the measurement. If the flashing symbol appears in three consecutive measurements, be sure to discuss this with your doctor.

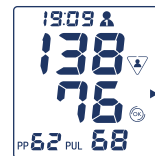
Symbol is not flashing: cardiac arrhythmias or irregular pulse rate detected.

The device is not a diagnostic instrument. It only indicates the potential presence of blood pressure-related conditions. **If this display should appear, it is essential to discuss this with your doctor to receive an accurate diagnosis.**

Caution: Make sure that you do not confuse the irregular pulse rate display  with the flashing pulse display , which shows the heart rate during each measurement (see page 51).

14. Motion control

Movements during the measurement may result in incorrect measurement values. The motion control registers stronger movements and these are indicated by the  symbol in the display. Other disruptive factors such as speaking, coughing or cuff noises can cause the symbol to be displayed. If this symbol appears, please repeat the measurement and avoid the disruptive factors.

The symbol is stored with the respective measurement.

15. Using the memory

The device has two memories (one per user) for storing measured values. Each memory can store up to 120 readings and the average value. The

measured results are automatically saved in the memory. If more than 120 measured values are saved, the oldest measured value (no.120) is deleted so that the latest value (no. 1) can be recorded.

Data retrieval

To retrieve data, press the Memory button of the user you wish to retrieve the data for while the device is switched off.

If the most recent measurement performed was a single measurement, the average value from the last three single measurements is calculated and displayed with "A" (Figure 1).

If the most recent measurement performed was a triple measurement, the average value of the last triple measurement is displayed with "A3" (Figure 2).

Press the Memory button again to display the most recently measured values (Figure 3). If the date/time function is activated, the memory number (Figure 3), the date (Figure 4) and the time (Figure 5) are displayed alternately.

Figure 3

Figure 4

Figure 5

Press the Memory button repeatedly to recall further measured values.


You can switch between the two memories by pressing Memory button 1 or 2.

The stored data remain visible in the display for around 30 seconds. The device then switches itself off.


Deleting data

To delete individual measured values, recall the measured value you wish to delete.

Press the Memory button again and hold down until the display firstly flashes and is then replaced with dashes.

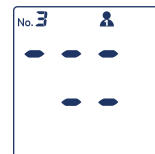To delete the entire memory, retrieve the respective average value (A or A3), press the Memory button again and hold down until the average value firstly flashes and is then replaced with dashes.

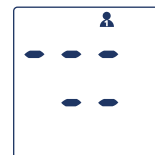16. Transferring the data to a PC

The device is equipped with a USB interface. With the help of the enclosed USB connection cable and the visomat double comfort software, you can transfer the measurement values from the blood pressure measurement device to your PC and evaluate them there.

The software has been specifically adapted to this blood pressure measurement device and can only be used with this device.

To operate the software, you need a PC with the Windows® operating system and a USB connection. You can find the current version of the visomat double comfort software as well as detailed instructions for installation and use at:

www.visomat.de/service/download

D What you should know about blood pressure

1. The systolic and diastolic blood pressure value

The circulatory system has the important task of supplying all organs and tissues in the body with sufficient blood, and of removing metabolites. The heart does this by contracting and expanding at a regular rhythm of around 60-80 times per minute. The pressure created by the blood flow on the artery walls when your heart contracts is referred to as systole. The pressure in the subsequent relaxation phase when the heart fills with blood again is referred to as diastole. In your daily measurement, you will measure both values.

2. Why you measure different values

Our blood pressure responds to internal and external influences like a sensitive measuring instrument. Affected by mental, physical and environmental influences, it varies all the time and never remains constant.

Reasons for fluctuating blood pressure values can be: moving, speaking, eating, consuming alcohol or nicotine, nervousness, inner tension, emotions, room temperature, recent urination or bowel movements, environmental influences such as movements and sounds, etc. Even changes in the weather and climate can have an effect on your blood pressure.

This explains why values measured at the doctor are often higher than those you obtain at home in your usual environment.

3. Why should you measure blood pressure regularly?

Even the time of day has an effect on your blood pressure. The values are at their highest during the day. In the course of the afternoon and in the evening, they drop slightly. They are low while you are sleeping, but rise again relatively quickly once you get up.

One-off and irregular measurements therefore say little about your actual blood pressure. A reliable assessment is only possible when measurements are taken regularly. Discuss the measured values with your doctor.

Technical information

1. Error messages

Error	Possible cause	Remedy
Display Err - 300	Excessive cuff pressure. The cuff was inflated to its maximum due to the movement of the arm or body.	-Repeat the measurement -Do not move your arm -Do not speak
	The air hose may be kinked or the air plug blocked.	Inspect the hose line and air plug for blockages
Display Err - 1	Measurement error; measurement could not be carried out due to movement or speaking.	-Repeat the measurement -Do not move your arm -Do not speak
	The acquired values lie outside the measuring range of the device (see "Technical data" on page 61).	Repeat the measurement
Display Err - 2	Inflation error; pressure built up too quickly or too slowly.	Check that the air connector is fitted properly, repeat measurement.
	Fault has been identified. The cuff is not correctly attached	-Repeat the measurement -Do not move your arm -Do not speak
Display Err - 3	Deflation rate too high or too low.	Call the Service hotline if this occurs repeatedly.
Display 0 0	Start/Stop button was accidentally pressed while the batteries were being replaced or while the mains adapter was being plugged in.	Switch unit off and then on again using the Start/Stop button.
 Symbol is not flashing	Arrhythmia, movement disorders, trembling, wobbling (objects), deep breathing etc.	Repeat measurement after resting for 3-5 minutes. Listen out for the device's pulse beep signals. Possible influences on the measurement results must be discussed with your doctor!
 Symbol is flashing	Atrial fibrillation	
	The motion control detected movement during the measurement.	-Repeat the measurement -Do not move your arm -Do not speak

Error	Possible cause	Remedy
Measured values/pulse values are extremely high or low.	Too little rest before taking the measurement. Movement or talking during the measurement.	Repeat the measurement after a break of 3-5 minutes. Pay attention to the usage information on page 45.
	Incorrect measurement posture	Repeat the measurement and maintain the correct posture. (Page 49)
Measurements taken in direct succession produce different measured values.	Blood pressure is not a fixed value. It can go up or down by more than 20 mmHg within a few minutes.	In order to obtain comparable results, always carry out the measurements under the same conditions.
Measurement was interrupted, the cuff deflated to 0 and the measurement restarted.	Device detects a fault or that the diastolic value is too low.	Device shows no fault. Repeat measurement
	Movement during measurement	If necessary, interrupt measurement and repeat after 5 minutes' rest.
No display after device is switched on	Batteries inserted incorrectly	Check position of batteries.
	Batteries flat	Replace batteries.
	Battery contacts dirty.	Clean battery contacts with a dry cloth.
Dashes in display after pressing Memory button	No measured values have been stored.	Repeat measurement
Measurement is interrupted	Batteries flat	Replace batteries.
The time is visible in the display, but the measurement is interrupted and the  symbol appears.	Batteries flat	Replace batteries.
Date and time are not visible in the display	Function was unintentionally deactivated	Activate the date and time function as described on page 46.

2. Customer service

Device repairs may only be performed by the manufacturer or an expressly authorised body. Please contact:

UEBE Medical GmbH
Bgm.-Kuhn-Str. 22
97900 Kilsheim, Germany
info@uebe.com
www.uebe.com

3. Technical data

Model:	visomat double comfort
Size:	L = 127.3 mm x W = 162.3 mm x H = 96 mm
Weight:	approx. 332 g without batteries and cuff
Display:	LCD display (liquid crystal display) for measured values and check displays
Memory:	2 memories each capable of storing 120 readings
Measurement method:	Combines auscultatoric and oscillometric determination of systole, diastole and pulse
Reference method of the clinical test:	Auscultatory measurement
Pressure display range:	0-300 mmHg
Measurement range:	Systolic: 50-250 mmHg Diastolic: 40-150 mmHg Pulse measurement: 30-199 pulses/minute
Measurement precision:	Blood pressure measurement (cuff pressure): ± 3 mmHg, Pulse rate: ± 5%
Power supply:	4 x 1.5 V AA LR6 alkaline batteries, durability: over 800 measurements in 2 years
	Optional: visomat mains adapter, output 6 V DC, 500 mA
Power consumption:	max. 4 W
Cuff:	Cuff type UWK for an arm circumference of 23-43 cm, applied part type BF

Operating conditions:	Ambient temperature 10°C to 40°C, relative humidity 15% to 85% non-condensing, air pressure 700 to 1060 hPa
Conditions for storage and transport:	Ambient temperature -20°C to +60°C Relative humidity 10% to 95%, non-condensing
Time to adjust from minimum/maximum storage temperature to operating temperature:	2 hours
Inflation and deflation:	Regulated automatically
IP rating:	IP 20: Protected against solid foreign bodies with a diameter from 12.5 mm, no protection against water.
Protection against electric shock:	Internal power supply, applied part type BF (cuff)
Expected service life:	5 years
Operating mode:	Continuous operation
Classification:	Internal power supply using batteries.

4. Original spare parts and accessories

The following original spare parts or accessories are available through specialist retailers:

- Universal cuff 23-43 cm type UWK
Part. no. 2405001, PZN-07514096
- Cuff 14-23 cm type USK
Part. no. 2405005, PZN-07514104
- visomat mains adapter
Part no. 2400020, PZN-11870684
- USB cable, Type A USB / Type Micro USB
Art. Nr. S2405690

- visomat double comfort software
download under www.visomat.de/service/download

Subject to technical modifications.

5. Metrological inspection

Generally, a metrological inspection is recommended every 2 years. However, professional users in Germany are obligated to do so in accordance with the Medical Devices EU Adaptation Act. This can be performed either by UEBE Medical GmbH, an authority responsible for metrology or an authorised maintenance service. Please refer to your national regulations.

Upon request, responsible authorities or authorised maintenance services receive a "Test instruction for metrological inspection" from the manufacturer.

Please only submit the device for metrological inspection together with the cuff and the instructions for use.

Important: No modifications, e.g. opening the device (except to replace the batteries), may be made to this device without the manufacturer's permission.

6. Explanation of symbols

This product complies with Regulation (EU) 2017/745 of the European Parliament and Council dated 5 April 2017 regarding medical products. It bears the marking CE 0123 (TÜV SÜD Product Service GmbH).

Consult the instructions for use

Caution

Manufacturer

Medical Device

Reference number = item number

Unique Device Identification:

(01): Device Identifier

(21): Production Identifier

Structure of the production identifier: YYWWMMMxxxxx

YY = Production year

WW = Production week (calendar week)

MMM = Model number

xxxxx = Serial number of individual device

Lot number

Temperature limitation

Humidity limitation

Keep dry

Degree of protection against electric shock TYPE BF

Protection class II (double insulation)

IP20

Protected against solid foreign bodies with a diameter from 12.5 mm, no protection against water.

USB connection socket

Direct current

Polarity inside plus

For indoor use only

7. Disposal

Technical equipment and batteries do not belong in household waste. They must be disposed of at appropriate collection and disposal points.

- Keep the device clean. Check for cleanliness after use. Please use a soft, dry cloth for cleaning. Do not use benzene, thinners or other strong solvents.
- The cuff may absorb sweat and other fluids. Inspect the cuff for marks or discolourations after each use. Please use a synthetic detergent for cleaning and carefully wipe the surface. Do not scrub or machine-wash. Air-dry thoroughly.
- Make sure that liquid cannot get into the air hose.
- Take care when storing the device. Ensure that no heavy objects are resting on the device or the cuff and that the air hose is not kinked. Do not wind the air hose too tightly.
- Do not apply force when folding or bending the cuff.
- To separate the cuff from the main device, please do not pull the air hose, but grip the air connector and gently pull it off.
- If the device is not to be operated for an extended period of time, please remove the batteries. Store the batteries out of the reach of children.

The device has been manufactured and tested with all due care. Nevertheless, in the event of defects upon delivery, we provide a warranty under the following conditions:

During the warranty period of 5 years from the date of purchase, we remedy such defects at our discretion and at our expense in our plant through repair work or replacement delivery of a defect-free device.

The warranty does not cover normal wear and tear of wearing parts or damage caused by failure to observe the instructions for use, improper handling (e.g. unsuitable power sources, breakage, leaking batteries) and/or disassembly of the device by the buyer. Furthermore, the warranty does not constitute grounds for asserting claims for damages against us.

Warranty claims can only be asserted during the warranty period and upon presentation of the proof of purchase. In the case of a warranty claim, the device along with the proof of purchase and a description of the complaint must be sent to:

UEBE Medical GmbH
Service-Center
Zum Läger 1
97900 Kulsheim, Germany

The cost of returning the device to our factory shall be borne by the sender. Complaints that are sent back without prepayment will not be accepted by UEBE.

The statutory claims and rights of the buyer against the seller (for example, defect claims, producer liability) are not restricted by this warranty.

Please note: In the case of a warranty claim, please make sure to enclose the proof of purchase.

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© Copyright 2025 UEBE Medical GmbH